

WakoVAC Vero


Vero细胞用细胞培养基

WakoVAC系列是本公司以GMP标准开发的用于疫苗生产的独特的细胞培养基。
目前正在开发的专门用于Vero细胞的培养基, 我们可以提供多个可用于测试的培养基原型。

■ 特点

- 不含动物源原料, 含有酵母水解物
- 广泛支持小试到大规模的疫苗生产
- 可以应对从研究到量产多种规模的使用
- 采用最新的GMP标准, 实施原材料变更管理

WakoVAC Vero的原型A以及原型B(图中培养基4、5), 按照以下培养条件, 和其他培养基相比具有最佳的细胞增殖效果。

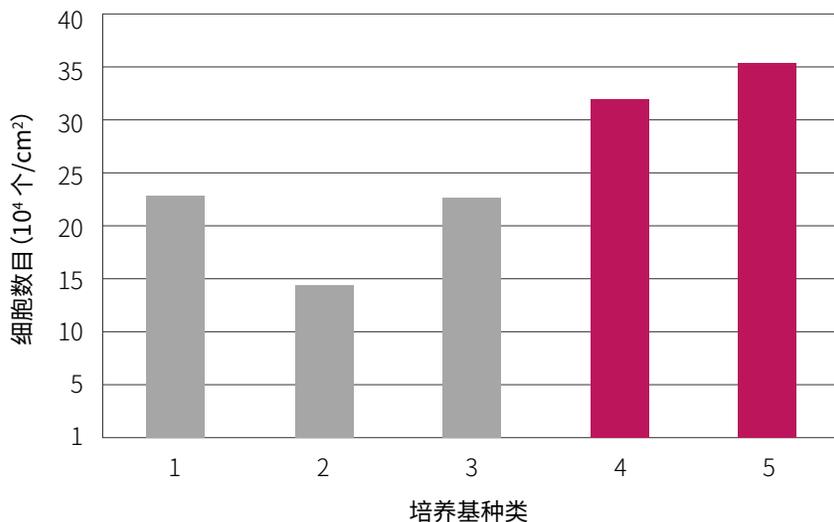

培养板: 12孔板 (康宁, #3513)

细胞类型: Vero细胞

细胞数: 30,000细胞/孔

培养条件: 37 °C、5%CO₂

培养时长: 5天

培养基用量: 2.0 mL

培养基:

1. EMEM/5% FBS
2. DMEM/5% FBS
3. 竞争品牌无血清培养基
4. WakoVAC Vero原型 A
5. WakoVAC Vero原型 B

本品可对企业用户提供试用装。

请联系富士胶片和光(广州)贸易有限公司在线客服或当地代理商, 进一步了解产品试用装申请事宜。

上述试剂仅供实验研究用, 不可用作“医药品”、“食品”、“临床诊断”等。

Listed products are intended for laboratory research use only, and not to be used for drug, food or human use. / Please visit our online catalog to search for other products from FUJIFILM Wako: <https://labchem-wako.fujifilm.com/> / This leaflet may contain products that cannot be exported to your country due to regulations. / Bulk quote requests for some products are welcomed. Please contact us.

富士胶片和光(广州)贸易有限公司

广州市越秀区先烈中路69号东山广场30楼3002-3003室

北京 Tel: 13611333218 上海 Tel: 021 62884751

广州 Tel: 020 87326381 香港 Tel: 852 27999019

询价: wkgz.info@fujifilm.com

官网: labchem.fujifilm-wako.com.cn

官方微信

目录价查询

